

anti- ARHGAP1 antibody

Product Information

Catalog No.: FNab00545

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

GTPase activator for the Rho, Rac and Cdc42 proteins, converting them to the putatively inactive GDP-bound state. Cdc42 seems to be the preferred substrate.

Immunogen information

Immunogen: Rho GTPase activating protein 1

Synonyms: CDC42GAP, RHOGAP1

Observed MW: 50 kDa UniprotID: Q07960

Application

Reactivity: Human

Tested Application: ELISA, WB, IHC

Recommended dilution: WB: 1:500-1:2000; IHC: 1:20-1:200

Image:


1

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 www.fn-test.com



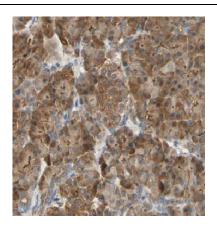

Immunohistochemistry of paraffin-embedded human pancreas tissue slide using FNab00545(ARHGAP1 Antibody) at dilution of 1:200

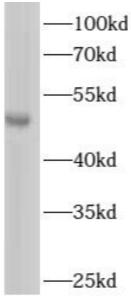

RAW 264.7 cells were subjected to SDS PAGE followed by western blot with FNab00545(ARHGAP1 Antibody) at dilution of 1:600