

anti- ATP5C1 antibody

Product Information

Catalog No.: FNab00704

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

Mitochondrial membrane ATP synthase(F(1)F(0) ATP synthase or Complex V) produces ATP from ADP in the presence of a proton gradient across the membrane which is generated by electron transport complexes of the respiratory chain. F-type ATPases consist of two structural domains, F(1)-containing the extramembraneous catalytic core, and F(0)-containing the membrane proton channel, linked together by a central stalk and a peripheral stalk. During catalysis, ATP synthesis in the catalytic domain of F(1) is coupled via a rotary mechanism of the central stalk subunits to proton translocation. Part of the complex F(1) domain and the central stalk which is part of the complex rotary element. The gamma subunit protrudes into the catalytic domain formed of alpha(3)beta(3). Rotation of the central stalk against the surrounding alpha(3)beta(3) subunits leads to hydrolysis of ATP in three separate catalytic sites on the beta subunits.

Immunogen information

Immunogen: ATP synthase, H+ transporting, mitochondrial F1 complex, gamma

polypeptide 1

Synonyms: ATP5C, ATP5C1, ATP5CL1, F ATPase gamma subunit

Observed MW: 35 kDa Uniprot ID: P36542

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>

1




Application


Reactivity: Human, Mouse, Rat Tested Application: ELISA, WB, IHC

Recommended dilution: WB: 1:500-1:2000; IHC: 1:100-1:500

Image:

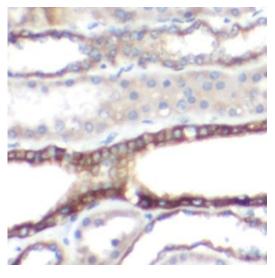

Immunohistochemistry of paraffin-embedded human kidney tissue slide using FNab00704(ATP5C1 Antibody) at dilution of 1:100

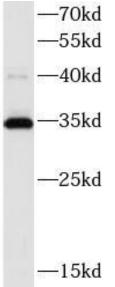

LO2 cells were subjected to SDS PAGE followed by western blot with FNab00704(ATP5C1 antibody) at dilution of 1:1001

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>