

anti- BBS4 antibody

Product Information

Catalog No.: FNab00818

Size: $100\mu g$ Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

The BBSome complex is thought to function as a coat complex required for sorting of specific membrane proteins to the primary cilia. The BBSome complex is required for ciliogenesis but is dispensable for centriolar satellite function. This ciliogenic function is mediated in part by the Rab8 GDP/GTP exchange factor, which localizes to the basal body and contacts the BBSome. Rab8(GTP) enters the primary cilium and promotes extension of the ciliary membrane. Firstly the BBSome associates with the ciliary membrane and binds to RAB3IP/Rabin8, the guanosyl exchange factor(GEF) for Rab8 and then the Rab8-GTP localizes to the cilium and promotes docking and fusion of carrier vesicles to the base of the ciliary membrane. The BBSome complex, together with the LTZL1, controls SMO ciliary trafficking and contributes to the sonic hedgehog(SHH) pathway regulation. Required for proper BBSome complex assembly and its ciliary localization. Required for microtubule anchoring at the centrosome but not for microtubule nucleation. May be required for the dynein-mediated transport of pericentriolar proteins to the centrosome.

Immunogen information

Immunogen: Bardet-Biedl syndrome 4

Synonyms: None

Observed MW: 50-60 kDa

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 www.fn-test.com

1




UniprotID: Q96RK4

Application


Reactivity: Human, Mouse, Rat Tested Application: ELISA, WB, IHC, IF

Recommended dilution: WB: 1:500-1:2000; IHC: 1:20-1:200; IF: 1:20-1:200

Image:

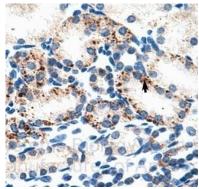

Immunohistochemistry of paraffin-embedded human kidney using FNab00818(BBS4 antibody) at dilution of 1:50

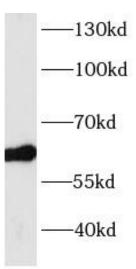

HEK-293 cells were subjected to SDS PAGE followed by western blot with FNab00818(BBS4 antibody) at dilution of 1:500