

anti- BRD7 antibody

Product Information

Catalog No.:	FNab00948
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	\geq 95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Acts both as coactivator and as corepressor.May play a role in chromatin remodeling.Activator of the Wnt signaling pathway in a DVL1-dependent manner by negatively regulating the GSK3B phosphotransferase activity.Induces dephosphorylation of GSK3B at 'Tyr-216'. Down-regulates TRIM24-mediated activation of transcriptional activation by AR(By similarity). Transcriptional corepressor that down-regulates the expression of target genes. Binds to target promoters, leading to increased histone H3 acetylation at 'Lys-9'(H3K9ac).Binds to the ESR1 promoter.Recruits BRCA1 and POU2F1 to the ESR1 promoter.Coactivator for TP53-mediated activation of target genes.Required for TP53-mediated cell-cycle arrest in response to oncogene activation. Promotes acetylation of TP53 at 'Lys-382', and thereby promotes efficient recruitment of TP53 to target promoters. Inhibits cell cycle progression from G1 to S phase.

Immunogen information

Immunogen:	bromodomain containing 7
Synonyms:	BP75, CELTIX1
Observed MW:	70-74kd
UniprotID :	Q9NPI1

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)


Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>


1



Application

Reactivity:Human, Mouse, RatTested Application:ELISA, WB, IHCRecommended dilution:WB: 1:500-1:2000; IHC: 1:20-1:200Image:

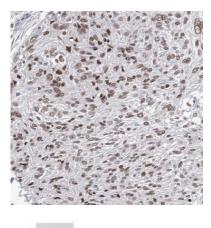

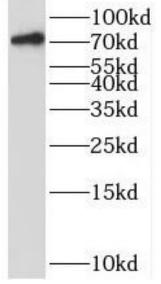

Immunohistochemistry of paraffin-embedded human cervical cancer tissue slide using FNab00948(BRD7 Antibody) at dilution of 1:50

HeLa cells were subjected to SDS PAGE followed by western blot with FNab00948(BRD7 antibody) at dilution of 1:500

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>