

anti- BRF2 antibody

Product Information

Catalog No.: FNab00955

Size: 100μg Form: liquid

Purification: Immunogen affinity purified

Purity: ≥95% as determined by SDS-PAGE

Host: Rabbit

Clonality: polyclonal

Clone ID: None IsoType: IgG

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

General activator of RNA polymerase III transcription. Factor exclusively required for RNA polymerase III transcription of genes with promoter elements upstream of the initiation sites. It recognizes a 50-60kd band in WB analysis.

Immunogen information

Immunogen: BRF2, subunit of RNA polymerase III transcription initiation factor,

BRF1-like

Synonyms: B related factor 2, BRF 2, BRF2, BRFU, FLJ11052, hBRFU, hTFIIIB50,

TFIIIB50

Observed MW: 50-60 kDa Uniprot ID: Q9HAW0

Application

Reactivity: Human, Mouse, Rat

Tested Application: ELISA, WB

Recommended dilution: WB: 1:500-1:1000

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>

1



Image:

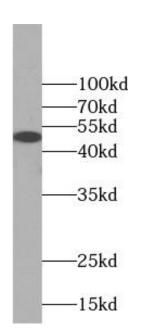

HT-1080 cells were subjected to SDS PAGE followed by western blot with FNab00955(BRF2 antibody) at dilution of 1:500

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel:(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>