

anti- CHRNB1 antibody

Product Information

Catalog No.:	FNab01682
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

After binding acetylcholine, the AChR responds by an extensive change in conformation that affects all subunits and leads to opening of an ion-conducting channel across the plasma membrane.

Immunogen information

Immunogen:	cholinergic receptor, nicotinic, beta 1(muscle)
Synonyms:	ACHRB, CHRNB
Observed MW:	50-57 kDa
UniprotID :	P11230

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

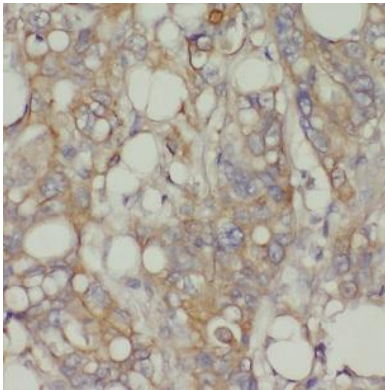

Immunohistochemistry of paraffin-embedded human colon cancer using FNab01682(CHRNB1 antibody) at dilution of 1:50

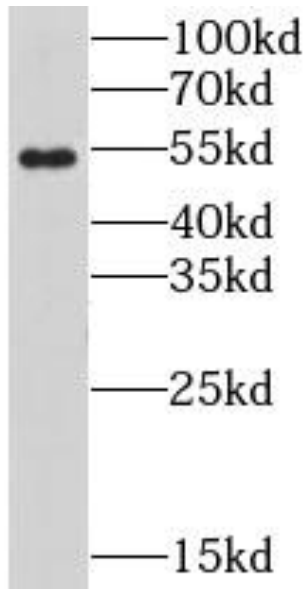

human liver tissue were subjected to SDS PAGE followed by western blot with FNab01682(CHRNB1 antibody) at dilution of 1:300