

anti- GGPS1 antibody

Product Information

Catalog No.:	FNab03445
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3 , -20°C for 24 months (Avoid repeated freeze / thaw cycles.)

Background

Catalyzes the trans-addition of the three molecules of IPP onto DMAPP to form geranylgeranyl pyrophosphate, an important precursor of carotenoids and geranylated proteins.

Immunogen information

Immunogen:	geranylgeranyl diphosphate synthase 1
Synonyms:	GGPP synthase, GGPPS, GGPPS1, GGPPSase, GGPS1
Observed MW:	35 kDa
Uniprot ID :	O95749

Application

Reactivity:	Human, Mouse, Rat; other species are not tested. Please decide the specificity by homology
Tested Application:	ELISA, WB, IHC, IP
Recommended dilution:	WB: 1:500-1:3000; IP: 1:500-1:1000; IHC: 1:20-1:200
Image:	


Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)


Tel : (0086)027-87384275

Fax: (0086)027-87800889

www.fn-test.com

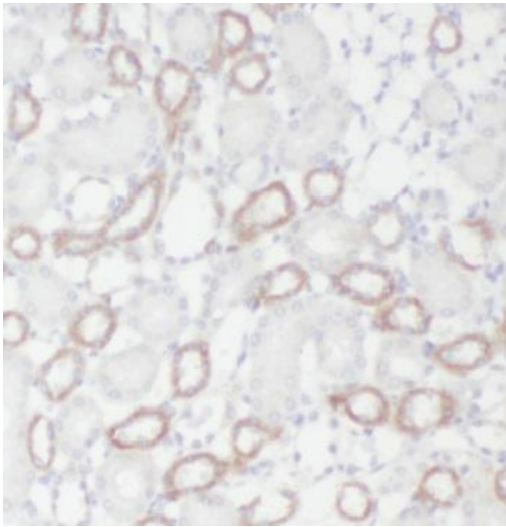

Immunohistochemistry of paraffin-embedded human kidney using FNab03445(GGPS1 antibody) at dilution of 1:100

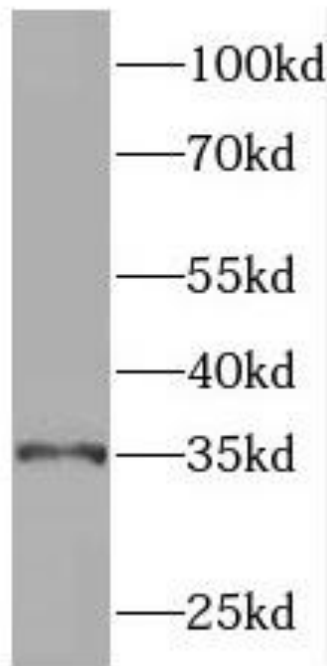

mouse brain tissue were subjected to SDS PAGE followed by western blot with FNab03445(GGPS1 antibody) at dilution of 1:500