

anti- GSDMD antibody

Product Information

Catalog No.:	FNab03670
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

GSDMD is a member of the gasdermin family. Members of this family appear to play a role in regulation of epithelial proliferation. GSDMD has been suggested to act as a tumor suppressor. Recently, GSDMD is identified as a new component of inflammasomes, and it is an executor of pyroptosis and required for interleukin-1 β secretion. Cleavage of GSDMD by inflammatory caspases determines pyroptotic cell death.

Immunogen information

Immunogen:	gasdermin D
Synonyms:	DF5L, DFNA5L, FKSG10, FLJ12150, gasdermin D, GSDMD, GSDMDC1
Observed MW:	50-53 kDa
Uniprot ID :	P57764

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, WB, IHC, IF

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Recommended dilution: WB: 1:200-1:2000; IHC: 1:20-1:200; IF: 1:10-1:100

Image:

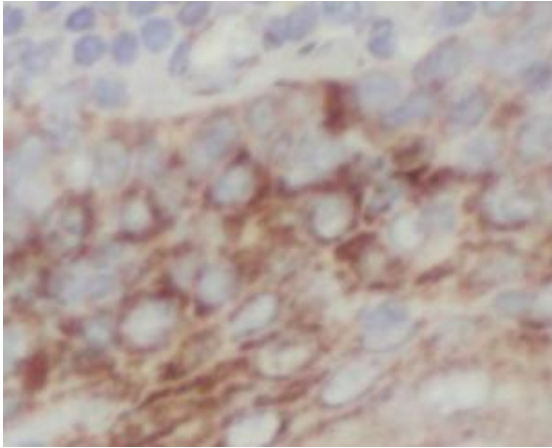

Immunohistochemistry of paraffin-embedded human skin using FNab03670(GSDMD antibody) at dilution of 1:50

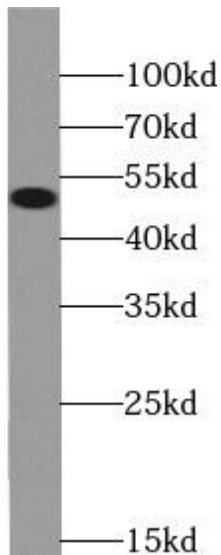

Jurkat cells were subjected to SDS PAGE followed by western blot with FNab03670(GSDMD Antibody) at dilution of 1:500

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889

www.fn-test.com