

## anti- HCLS1 antibody

### Product Information

Catalog No.:	FNab03786
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

### Background

HCLS1 (hematopoietic cell-specific Lyn substrate 1) is a protein-coding gene. Diseases associated with HCLS1 include chronic lymphocytic leukemia. GO annotations related to this gene include protein kinase binding and protein complex binding. An important paralog of this gene is DBN1.

### Immunogen information

Immunogen:	hematopoietic cell-specific Lyn substrate 1
Synonyms:	HS1
Observed MW:	72 kDa
UniprotID :	P14317

### Application

Reactivity:	Human, Mouse
Tested Application:	ELISA, WB, IHC
Recommended dilution:	WB: 1:500 - 1:2000; IHC: 1:50 - 1:200


### Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)


Tel : (0086)027-87384275

Fax: (0086)027-87800889 [www.fn-test.com](http://www.fn-test.com)

Image:

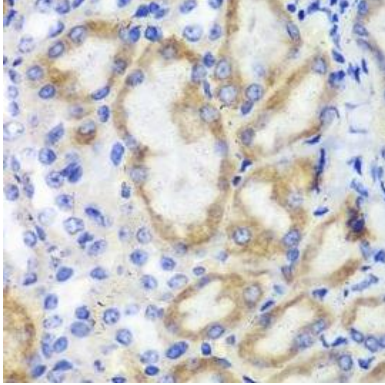

Immunohistochemistry of paraffin-embedded mouse kidney using FNab03786( HCLS1 Antibody) at dilution of 1:50

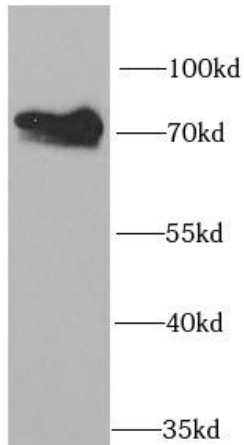

Jurkat cells were subjected to SDS PAGE followed by western blot with FNab03786( HCLS1 Antibody) at dilution of 1:1000

**Wuhan Fine Biotech Co., Ltd.**

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 [www.fn-test.com](http://www.fn-test.com)