

anti- HNRNPF antibody

Product Information

| | |
|---------------|---|
| Catalog No.: | FNab03952 |
| Size: | 100µg |
| Form: | liquid |
| Purification: | Immunogen affinity purified |
| Purity: | ≥95% as determined by SDS-PAGE |
| Host: | Rabbit |
| Clonality: | polyclonal |
| Clone ID: | None |
| IsoType: | IgG |
| Storage: | PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.) |

Background

Component of the heterogeneous nuclear ribonucleoprotein(hnRNP) complexes which provide the substrate for the processing events that pre-mRNAs undergo before becoming functional, translatable mRNAs in the cytoplasm. Plays a role in the regulation of alternative splicing events. Binds G-rich sequences in pre-mRNAs and keeps target RNA in an unfolded state.

Immunogen information

| | |
|--------------|---|
| Immunogen: | heterogeneous nuclear ribonucleoprotein F |
| Synonyms: | HNRPF |
| Observed MW: | 46 kDa |
| UniprotID : | P52597 |

Application

| | |
|--------------------------|-------------------------------|
| Reactivity: | Human, Mouse, Rat |
| Tested Application: | ELISA, WB, IHC |
| Recommended dilution:WB: | 1:500-1:2000; IHC: 1:20-1:200 |


Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)


Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

Image:

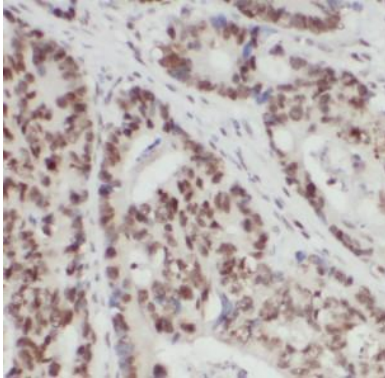

Immunohistochemistry of paraffin-embedded human colon cancer using FNab03952(HNRNPF antibody) at dilution of 1:50

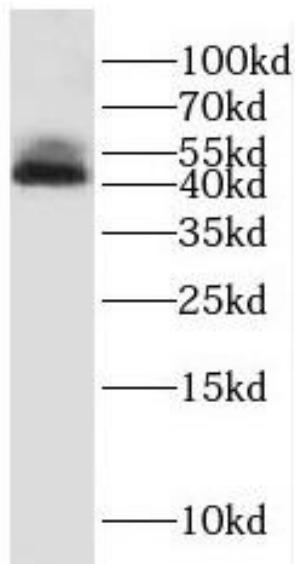

HeLa cells were subjected to SDS PAGE followed by western blot with FNab03952(HNRNPF antibody) at dilution of 1:1000

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve. East Lake High-Tech Development Zone. Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com