

anti- LONP1 antibody

Product Information

Catalog No.: FNab04821

Size: 100μg Form: liquid

Purification: Protein A+G purification

Purity: ≥95% as determined by SDS-PAGE

Host: Mouse

Clonality: monoclonal

Clone ID: 9H10 IsoType: IgG2a

Storage: PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12

months (Avoid repeated freeze / thaw cycles.)

Background

ATP-dependent serine protease that mediates the selective degradation of misfolded, unassembled or oxidatively damaged polypeptides as well as certain short-lived regulatory proteins in the mitochondrial matrix. May also have a chaperone function in the assembly of inner membrane protein complexes. Participates in the regulation of mitochondrial gene expression and in the maintenance of the integrity of the mitochondrial genome. Binds to mitochondrial promoters and RNA in a single-stranded, site-specific, and strand-specific manner. May regulate mitochondrial DNA replication and/or gene expression using site-specific, single-stranded DNA binding to target the degradation of regulatory proteins binding to adjacent sites in mitochondrial promoters. Endogenous substrates include mitochondrial steroidogenic acute regulatory(StAR) protein.

Immunogen information

Immunogen: lon peptidase 1, mitochondrial

Synonyms: PRSS15
Observed MW: 100 kDa
Uniprot ID: P36776

1

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave. East Lake High-Tech Development Zone. Wuhan, Hubei, China (430206)

Tel :(0086)027-87384275 Fax: (0086)027-87800889 <u>www.fn-test.com</u>




Application


Reactivity: Human, Mouse, Rat Tested Application: ELISA, WB, IHC, IF

Recommended dilution: WB: 1:500-1:10000; IHC: 1:400-1:1600; IF: 1:10-1:100

Image:

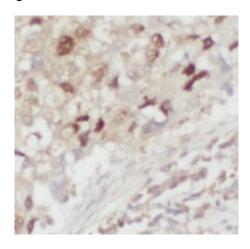

Immunohistochemistry of paraffin-embedded human lung cancer tissue slide using FNab04821(LONP1 Antibody) at dilution of 1:800 heat mediated antigen retrieved with Tris-EDTA buffer(pH9).

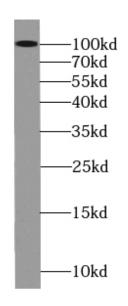

Jurkat cells were subjected to SDS PAGE followed by western blot with FNab04821(LONP1 Antibody) at dilution of 1:5000