

anti- LRP2BP antibody

Product Information

Catalog No.:	FNab04846
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

May act as an adapter that regulates LRP2 function.

Immunogen information

Immunogen:	LRP2 binding protein
Synonyms:	KIAA1325
Observed MW:	45 kDa
UniprotID :	Q9P2M1

Application


Reactivity:	Human, Rat
Tested Application:	ELISA, WB
Recommended dilution:WB:	1:500-1:2000
Image:	

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

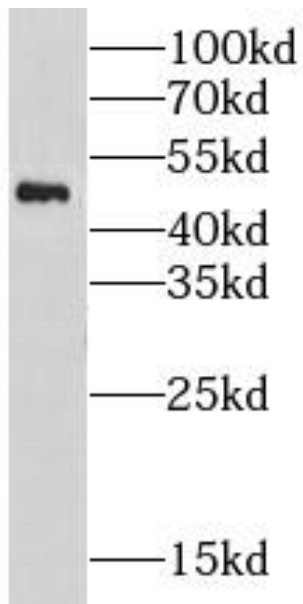

A549 cells were subjected to SDS PAGE followed by western blot with FNab04846(LRP2BP antibody) at dilution of 1:300