

anti- LRRC41 antibody

Product Information

Catalog No.:	FNab04857
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

Probable substrate recognition component of an ECS(Elongin BC-CUL2/5-SOCS-box protein) E3 ubiquitin ligase complex which mediates the ubiquitination and subsequent proteasomal degradation of target proteins.

Immunogen information

Immunogen:	leucine rich repeat containing 41
Synonyms:	MUF1
Observed MW:	80-90 kDa
UniprotID :	Q15345

Application


Reactivity:	Human, Mouse, Rat
Tested Application:	ELISA, IHC, WB
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:20-1:200
Image:	

Wuhan Fine Biotech Co., Ltd.


B9 Bld, High-Tech Medical Devices Park, No. 818 GaoxinAve.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel : (0086)027-87384275

Fax: (0086)027-87800889 www.fn-test.com

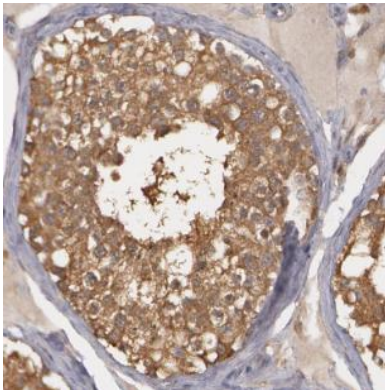

Immunohistochemistry of paraffin-embedded human testis using FNab04857(LRRC41 antibody) at dilution of 1:50

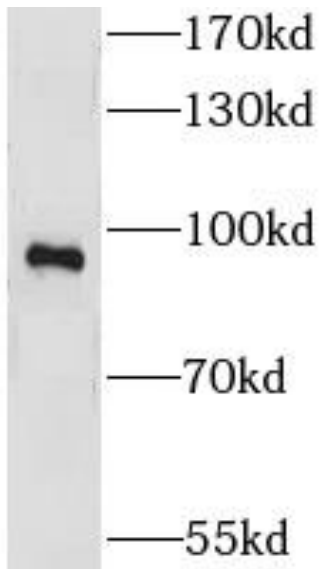

human brain tissue were subjected to SDS PAGE followed by western blot with FNab04857(LRRC41 antibody) at dilution of 1:100