

anti- ATF4 antibody

Product Information

Catalog No.:	FNab10416
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	\geq 95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months (Avoid repeated freeze / thaw cycles.)

Background

The antibody recognizes the 38kd ATF4 protein and its phosphorylated forms (50kd). ATF4 can bind DNA as a homodimer and as a heterodimer.

Immunogen information

Immunogen:	activating transcription factor 4(tax-responsive enhancer element B67)
Synonyms:	ATF4, CREB 2, CREB2, DNA binding protein TAXREB67, TAXREB67, TXREB
Observed MW:	50-58 kDa, 35-38 kDa
Uniprot ID :	P18848


Application

Reactivity:	Human		
Tested Application:	ELISA, WB, IHC		
Recommended dilution: WB: 1:500 - 1:2000; IHC: 1:50 - 1:200			
Image:			


Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206) Tel :(0086)027-87384275



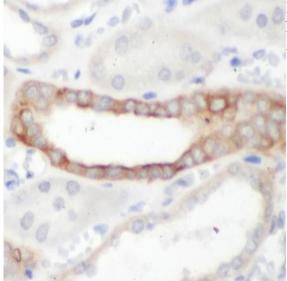

Immunohistochemistry of paraffin-embedded human kidney tissue slide using FNab10416(ATF4 Antibody) at dilution of 1:50

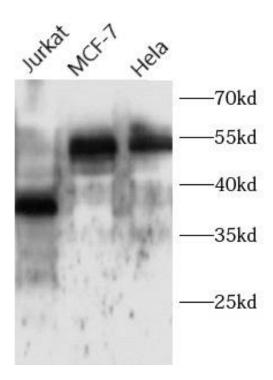

Various lysates were subjected to SDS PAGE followed by western blot with FNab10416(ATF4 antibody) at dilution of 1:1000

Wuhan Fine Biotech Co., Ltd.

B9 Bld, High-Tech Medical Devices Park, No. 818 Gaoxin Ave.East Lake High-Tech Development Zone.Wuhan, Hubei, China(430206)

Tel :(0086)027-87384275

www.fn-test.com